


Revolution Mobile Curtain

REV-PMC

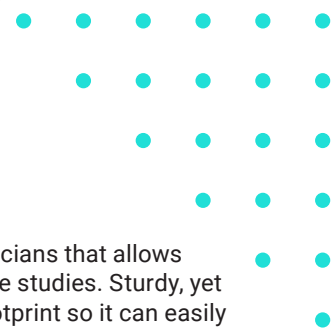

A unique new design for operating physicians that allows unparalleled protection for close-to-table studies. Sturdy, yet thoughtfully engineered with a small footprint so it can easily be positioned exactly where you want it for maximum scatter radiation protection.

- Special 4 caster T-bar base design to stabilize the unit
- Smooth easy maneuvering
- Special spring mechanism allows for easy raising, lowering, and automatically holding shield
- .50 mm LE double protective curtain
- 2 high-quality hospital grade locking casters

Model Number	Description
REV-PMC	Revolution Mobile Curtain
Dimensions	Adjustable height 43½" - 66" Base: 34" W x 25" D
Curtain	38"H x 30"W, Revolution Grey Custom colors available.
Frame Colors	White, Grey, or Black Custom colors available. Ships in standard white unless specified.
Frame Options	Powder Coated Steel Stainless Steel (Upgrade)
Wheels	4 swivel casters, 2 with brakes
Crate Size	Box 1: 13" x 42" x 9" Box 2: 26" x 30" x 8"
Crate Weight	Box 1: 45 lbs / 20.41 kg Box 2: 30 lbs / 13.60 kg

Spring mechanism makes it easy to raise and lower

